

Safe-T-Vue® 10 • Product No. 7201

Quality Assurance: Safe-T-Vue® 10 Non-Reversible Blood Bag Temperature Indicator Product No. 7201

PLEASE NOTE: ALL BOX OF 50 INFORMATION AND FOIL LID COPY FOR SAFE-T-VUE® 10 IS PRINTED IN BLACK.

Test performed on Random Sampling of Finished Product Lot # S10 011608

Expiration Date: 2010-01

Non - Chemistry

Mark ✓ for "Accepted" and X for "Not Accepted."

Box	Adh. Label	Foil Lid	Close Function	Fill Level
✓ Copy	✓ Copy	✓ Copy	✓ Formed	✓ Acceptable
✓ Tack	✓ Tack	✓ Position	✓ Locked/Stable	Underfill
✓ Clean	✓ Color	✓ Clean Peel	Hinge	Overfill
✓ Box Damage	Tray	Disk	✓ Elasticity	Color Label
	✓ Edges	✓ Appearance		✓ Color
	✓ Color	✓ Fit		

Test Protocol – Finished Goods Quality Assurance

Test performed on random sampling of finished product per ANSI ASQZ1.4-2003

Performed by using a "DigiSense" Thermocouple Electronic Thermometer. Using a cold (2° to 6° C) bag of simulated blood (a glycerol-water-red dye solution) with the electric thermometer placed and sealed into the core of the blood bag to permit observation of temperature changes. For detailed protocol request QAP 0300-01 REV 0 Mar 01, 2005 [Incoming inspection of 7201 Safe-T-Vue® 10]

Temperature Change and Corresponding Safe-T-Vue® Change

Sample Set 1	DATA	Sample Set 2
1 to 4 °C	<i>Starting Temperatures</i>	1 to 4 °C

OBSERVATIONS


Temperature	Condition
9.2 °C	Red Spots - Random WHITE CENTER
9.8 °C	Red Filling In WHITE CENTER
10.1 °C	Full Red Color RED CENTER

Temperature	Condition
9.2 °C	Red Spots – Random WHITE CENTER
9.8 °C	Red Filling In WHITE CENTER
10.1 °C	Full Red Color RED CENTER

Temperature accuracy is +/-0.3°C.

Results of Testing Lot #S10 011608; Test Date: January 25, 2008

Red Spots - Random 9.2°C; Red Filling in White Center 9.8°C; Full Red Color at 10.1°C.

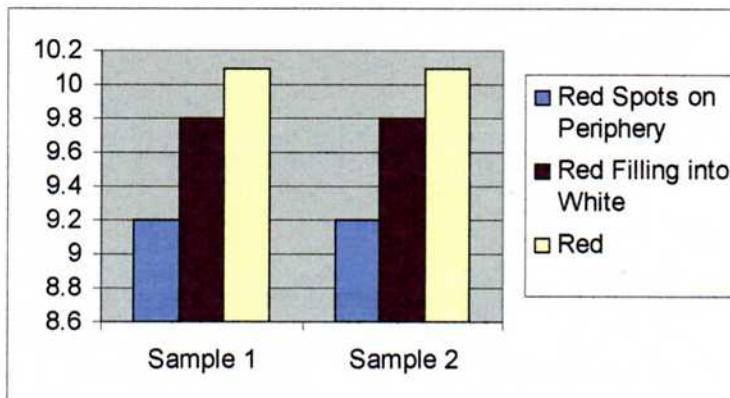

QUALITY ASSURANCE: APPROVED LOT # S10 011608 ; Edward N. Sharpless **January 25, 2008**

Laboratory Director

Date