

Safe-T-Vue® 10 • Product No. 7201

Quality Assurance: Safe-T-Vue® 10 Non-Reversible Blood Bag Temperature Indicator Product No. 7201
 PLEASE NOTE: ALL BOX OF 50 INFORMATION AND FOIL LID COPY FOR SAFE-T-VUE® 10 IS PRINTED IN BLACK.

Lot # S10 033011
 Expiration Date: 2013-03

Non - Chemistry

Mark ✓ for "Accepted" and X for "Not Accepted."

Box	Adh. Label	Foil Lid	Close Function	Fill Level
✓ Copy	✓ Copy	✓ Copy	✓ Formed	✓ Acceptable
✓ Tack	✓ Tack	✓ Position	✓ Locked/Stable	Underfill
✓ Clean	✓ Color	✓ Clean Peel	Hinge	Overfill
✓ Box Damage	Tray	Disk	✓ Elasticity	Color Label
	✓ Edges	✓ Appearance		✓ Color
	✓ Color	✓ Fit		

Test Protocol – Finished Goods Quality Assurance

Test performed on random sampling of finished product per ANSI ASQZ1.4-2003

Data from a "DigiSense" Thermister Electronic Thermometer, using a cold (start temperature 1° to 3° C) bag of simulated blood (10% glycerol-water) with thermometer probe positioned in center of blood bag. Agitation of blood bag necessary for collection of accurate temperature-color change data. Download detailed protocol, from www.williamlabs.com, QAP 0300-01 REV 0 "Receiving Inspection – Temperature Indicators".

Temperature Change and Corresponding Safe-T-Vue® Change

Sample Set 1	DATA	Sample Set 2
1 to 4 °C	<i>Starting Temperatures</i>	1 to 4 °C


OBSERVATIONS

Temperature	Condition
9.1 °C	Red Spots - Random WHITE CENTER
9.7 °C	Red Filling In WHITE CENTER
10.0 °C	Full Red Color RED CENTER

Temperature	Condition
9.1 °C	Red Spots – Random WHITE CENTER
9.6 °C	Red Filling In WHITE CENTER
10.0 °C	Full Red Color RED CENTER

Temperature accuracy is +/-0.3°C.

Results of Testing: Red Spots - Random 9.1°C; Red Filling in White Center 9.7°C; Full Red Color at 10.0°C.

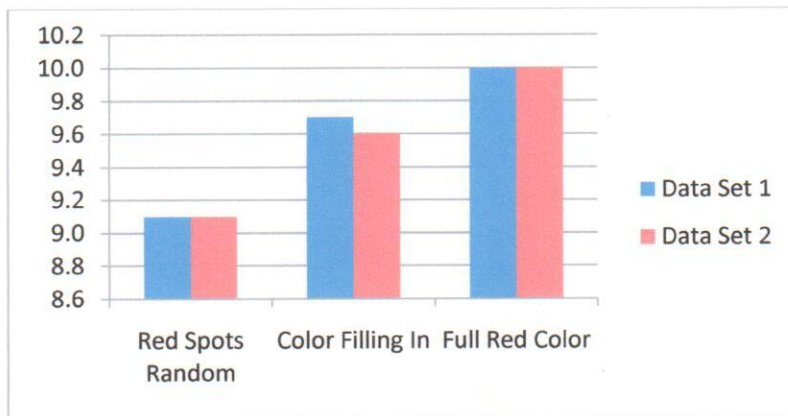

QUALITY ASSURANCE: APPROVED LOT # S10 033011 ; Edward N. Sharpless APR 27, 2011
 Laboratory Director Date